


5mm Conductivity Electrode

NEW Small Diameter EC electrode for MANTECH's MT Series Environmental Titration Systems

Part Number Information

PCE-96-CT1010

Fast Response Conductivity Cell, Compact 5mm Diameter, 2 Metre Cable, Glass, PT. Type with ATC, K=1, 25mm cell length, Quick Response. This new small diameter conductivity electrode is ideal for low sample volumes and skinny tubes. Takes less space in vials than standard 12mm conductivity electrodes. Well suited for performing conductivity measurements in 15mL vials, or combining with pH and Alkalinity directly in a standard 50mL sample tube. Thermistor probe built into electrode for automated temperature reading and compensation. Compatible with all MANTECH Analyzers currently measuring conductivity (note that a probe holder change will be required).

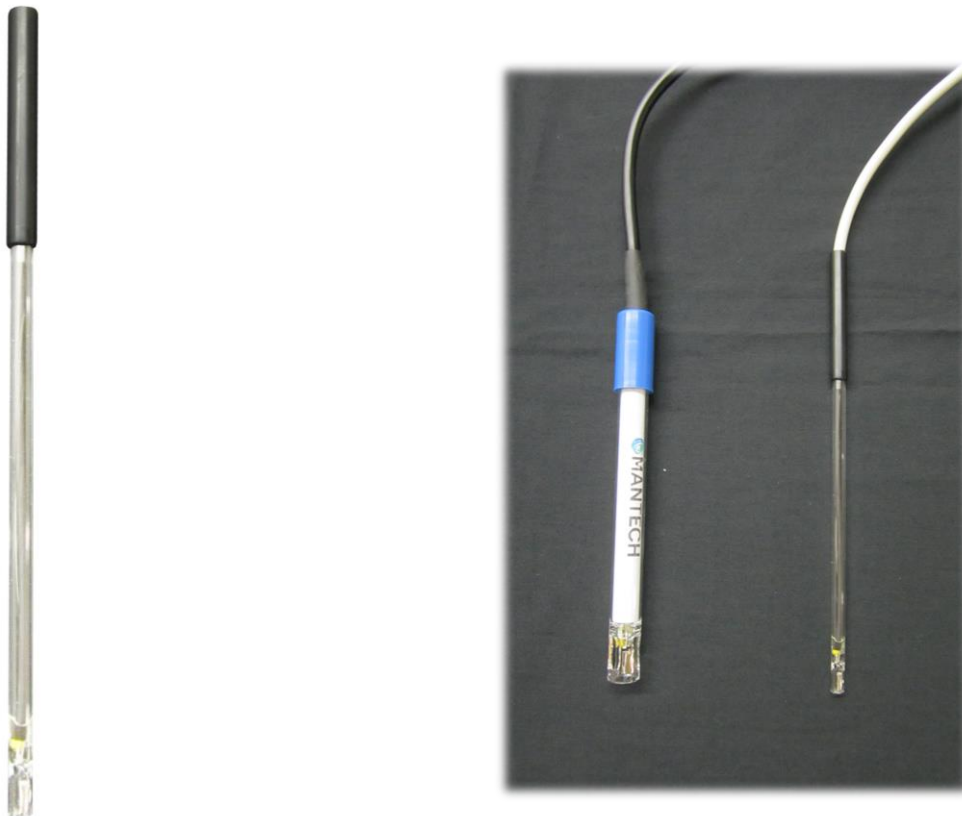

Pictured with 12mm EC probe on left, 5mm EC probe on right