



# MANTECH

OPTIMIZE YOUR RESULTS. PROTECT OUR ENVIRONMENT.



**5-DAY BIOCHEMICAL OXYGEN  
DEMAND (BOD5) SOFTWARE PACKAGE**



# COMPLIANCE MEETS FLEXIBILITY

BOD Pro™ is lab automation software coded and designed for MANTECH analyzers. Offered with automation or as a stand-alone software, BOD Pro™ is built with features to simplify and streamline 5-day biochemical oxygen demand (BOD5) analyses while delivering valid and compliant results.

**SOFTWARE ONLY PACKAGE**  
**\$2,090 USD**



- BOD PRO™ CONFORMS TO THE FOLLOWING REGULATORY NORMS:**
- ASTM STANDARD METHODS 5210B, 23<sup>RD</sup> AND 24<sup>TH</sup> EDITIONS
  - EPA 405.1
  - ISO EN 5815-1
  - ISO EN 1899-1, -2

## Probe Options

Membrane or optical

Single, dual or triple probes

YSI MultiLab 4010-W integration or manual entry required

## SOFTWARE FEATURES



### Customizable Shortcut Buttons

- Simplifies and reduces prep on routine sample runs, maintenance setup and more



### IoT Device Connectivity

- Bi-directional communication with hardware and software



### Import & Export LIMS Capabilities

- Import sample run information
- Auto-exporting feature to LIMS or network location



### Enable OS & Security Updates

- Designed for Windows 10 and 11 Pro OS

# UPGRADE TO BOD PRO™ TODAY!

**Plan on upgrading to Windows 10 or 11 Pro?** Due to Microsoft's major security updates, PC-BOD is not guaranteed to be operational with the latest OS updates. MANTECH's **NEWEST** BOD Pro™ software is designed for Windows 10 & 11 Pro OS for enhanced cyber security and dependability.



## OPTIMIZED SECURITY

- Flexible user account system with easy configuration and user defined permissions



## ENHANCED USER OPERATION

- Intuitively designed with colour coding, visual rack display and live script monitoring features



## HARDWARE INTEGRATION

- PC-BOD hardware is not compatible with BOD Pro™ software and requires MANTECH's [LATEST hardware](#)

## CONFIGURATIONS TO OPTIMIZE YOUR THROUGHPUT

### SIMPLIFY MANUAL DO READINGS




- ✓ Automated calculations
- ✓ Incubation period tracking
- ✓ Use your existing DO meter

### AUTOMATE DO READINGS AND SAMPLE PREPARATION




- ✓ Up to 108 position autosamplers
- ✓ Optional dilution water, seed and/or inhibitor addition


# BOD PRO™ SOFTWARE

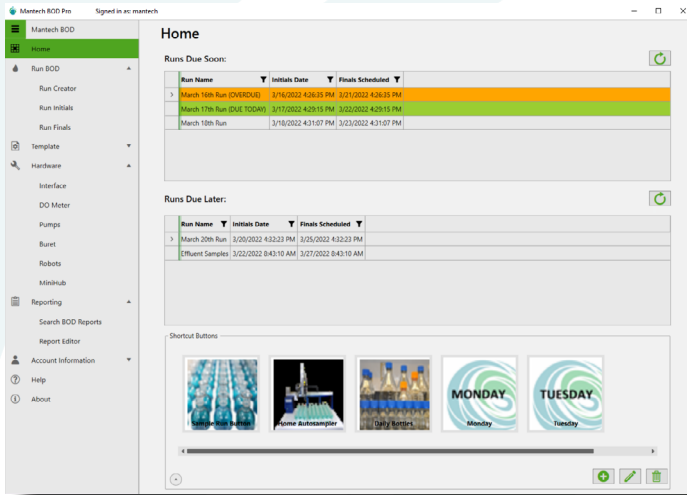


Home Screen With Customizable Shortcut Buttons

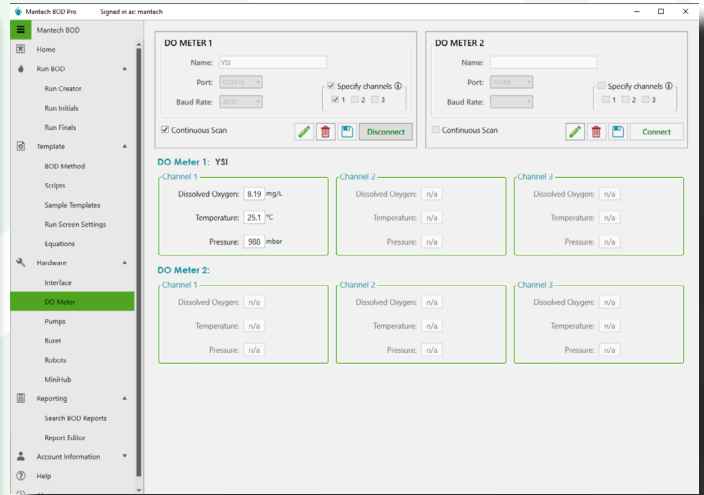

Up to Two Single, Dual or Triple Probe DO Meters

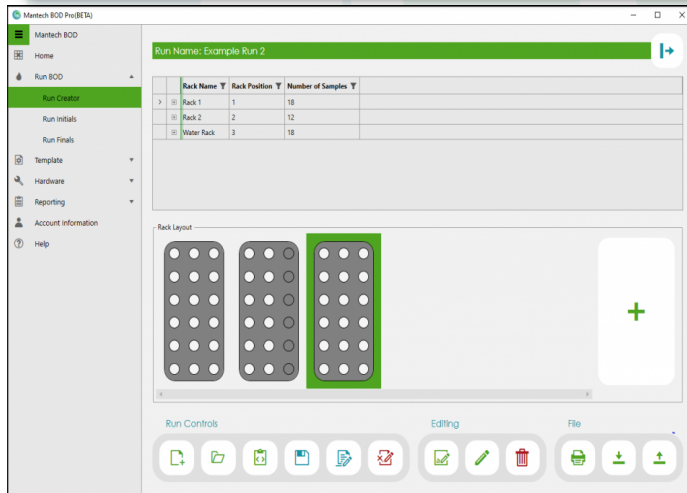

Build Run Templates With Intuitive Rack Display

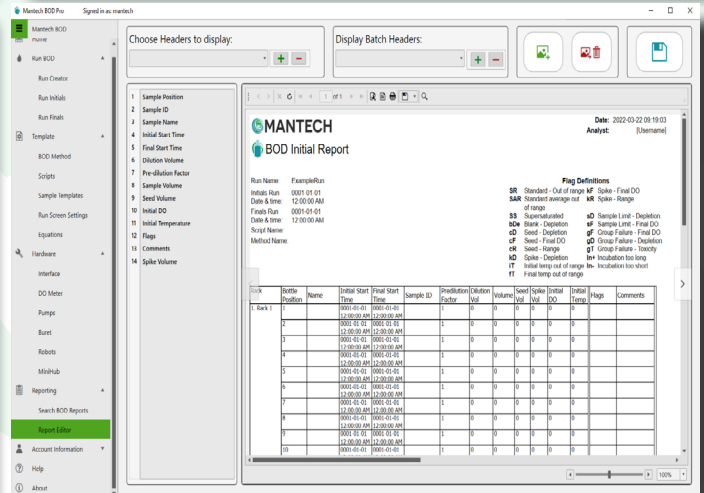

Add Up to 10 User-Defined Columns in Report Templates



5473 HIGHWAY 6 NORTH  
GUELPH, ONTARIO, CANADA  
N1H 6J2

WWW.MANTECH-INC.COM  
INFO@MANTECH-INC.COM  
+1 (519) 763-4245  
+1 (866) 763-2122

@mantech\_inc  
 linkedin.com/company/mantech-inc

© 2025 MANTECH. All rights reserved.  
PeCOD® is a registered trademark of  
MANTECH.